

UDS-P LED



Comfortable scaling

product details

A good helper for dentist

Optical handpiece, more convenient for clinical operation.

Automatic frequency tracking ensures that the machine always works on the best frequency and performs more steadily.

The handpiece is detachable and can be autoclaved to the high temperature of 135°C and the pressure of 0.22MPa.

Digitally controlled, easy operation and more efficient for scaling.




Detachable handpiece with light

Shell of the handpiece is made of special imported plastic from America.


Without cracking or deformation under the pressure.

Detachable handpiece with light, better visibility.


Note: This type can match with handpiece without light.

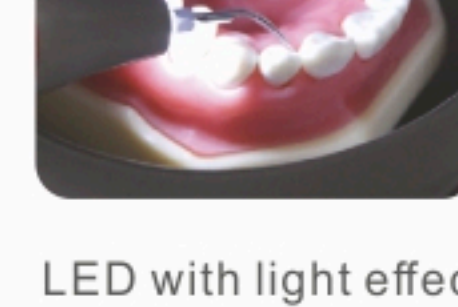

LED with light effect

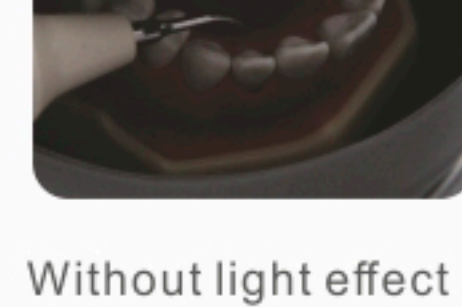

Without light effect

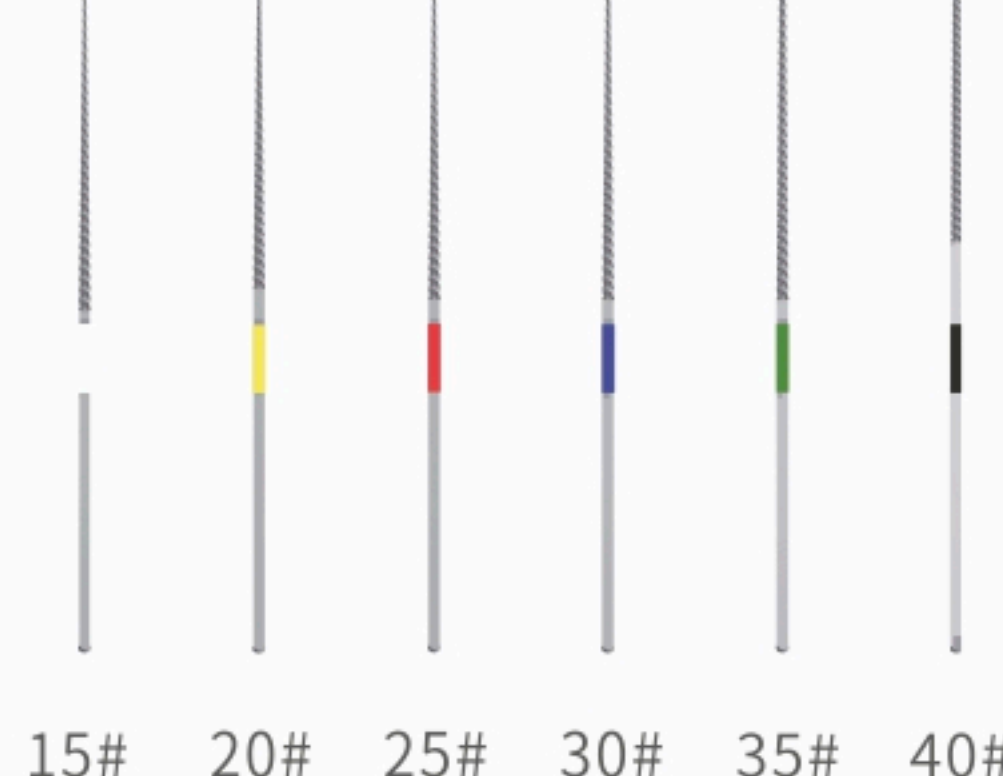

15# 20# 25# 30# 35# 40#

NITI-File (included)

15# 20# 25# 30#

TIPS

> Supragingival scaling



> Implant maintenance




P90
Implant maintenance
Clean implants and restorations

Subgingival scaling



Root canal treatment

Endo scaling

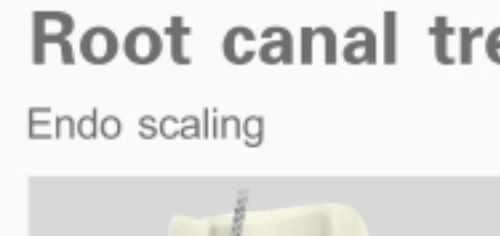


For removing the calcification and bad fillings in root canal




For removing the calcification in the coronal third of the root canal



Root canal retreatment, removing foreign body

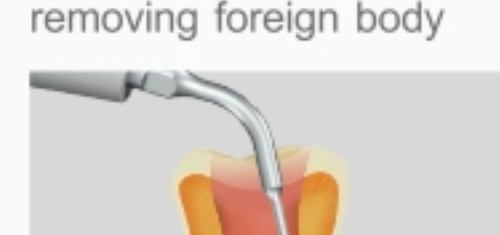

Looking for root canal orifice

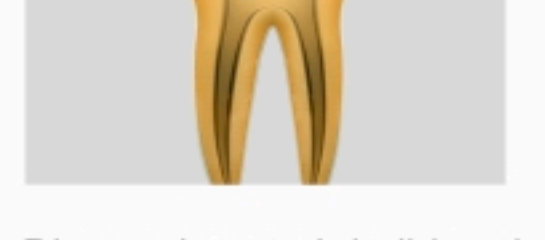

Special material [E6, E7], extremely thin and soft, good conductivity to ultrasound

For removing foreign body in 1/3 deep root canal

For removing foreign body in 1/3 root canal

For removing foreign body and broken instruments in root canal (stainless steel)

For removing root canal fillings and broken instruments (stainless steel)

Diamond-coated, ball head, for looking for root canal orifice with security, accuracy and efficiency

Cavity preparation



G30 For the preparation between pits and fissures

G31 For the preparation of mesioocclusion cavity

G32 For the preparation of distal mesial cavity

G30

G31

G32

SB1

SB2

SB3